

A VARIETY OF CAPSULORHEXIS FORCEPS


STANDARD

3.0 mm


4-03315T Utrata Capsulorhexis Forceps with Scale

Cystotome tips, TC*

Scale for measuring the size of capsulorhexis

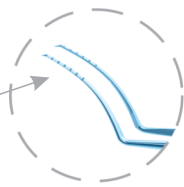

Curved shafts prevent corneal deformation

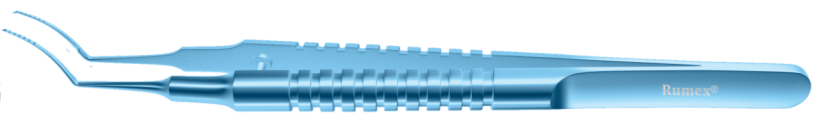

4-0331D Disposable Utrata Capsulorhexis Forceps

STERILE



Cystotome tips




- ✓ Meets quality standards applied to reusable tools
- ✓ Always ready to use



SMALL-INCISION

2.0 mm

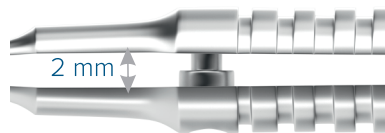


Internal limiter prevents the jaws from opening wider than intended

4-0312S Microcoaxial Capsulorhexis Forceps with Alignment Mechanism

Cystotome tips, TC*

Limited to 2.00 mm

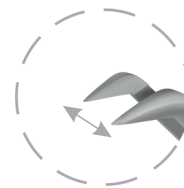

4-032S Small-Incision Capsulorhexis Forceps with Limiter and Wide Opening

Cystotome tips, TC*

Scale for measuring the size of capsulorhexis

Limited to 2.00 mm



View-port for better visualization



*TC – Tungsten Carbide coating makes the platform rough enough for controlled gripping and prevents adhesion of the capsule.

A VARIETY OF CAPSULORHEXIS FORCEPS

SMALL-INCISION


1.5 mm

4-0395

Capsulorhexis Forceps with Cross-Action and Scale

Cystotome tips, TC*

Scale for measuring the size of capsulorhexis

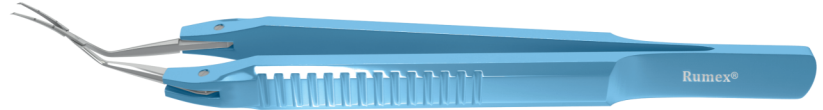

- ✓ Cross-action prevents hyperextension of the incision
- ✓ Hinge joints ensure even distribution of pressure at three points of forceps

4-033S

Small-Incision Capsulorhexis Forceps with Double Cross-Action and Scale

Cystotome tips, TC*

Scale for measuring the size of capsulorhexis



- ✓ Double cross-action guarantees fixed opening for free movements within the wound

MICS SUB

sub 1.0 mm

4-03742

Lesieur Capsulorhexis Forceps with Internal Ruler

12-003T


Universal Squeeze Handle

Scale for measuring the size of capsulorhexis



Anti-glare coating

Cystotome tips

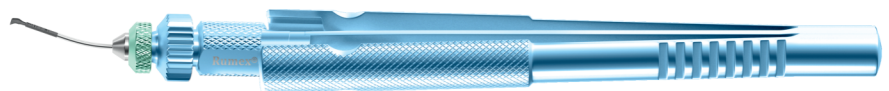

- ✓ 23 Ga tube minimally disrupts biomechanics of cornea
- ✓ MICS tips are used with RUMEX Universal Handle and can be easily disassembled for cleaning after the surgery

*TC – Tungsten Carbide coating makes the platform rough enough for controlled gripping and prevents adhesion of the capsule.

RUMEX INTERNATIONAL COMPANY

Head Office: 14240 Carlson Circle,
Building K, Suite 8, Tampa,
FL 33626, USA

 www.rumex.com

USA & CANADA

 +1 (727) 535 9600
+1 (877) 77 RUMEX (Toll-Free)
 +1 (727) 535 8300
 usacs@rumex.com

EUROPE, ASIA, AFRICA, LATIN AMERICA

 +371 6616 3182
 customerservice@rumex.com